




Wheelchair Accessible Vehicle **Securement**

How-To Guide

Don't rush!
Take your time.

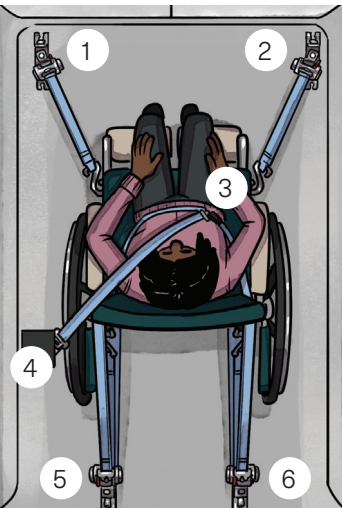

Use all six securement straps. Do not use broken equipment.



Explain to your passengers what you are doing.

I'm going to secure the shoulder belt to the lap belt...



Listen to your passengers.

Can you make this lap belt a little looser?




Keep your vehicle clean.


Onboarding




1 Slide front securement into tracks.

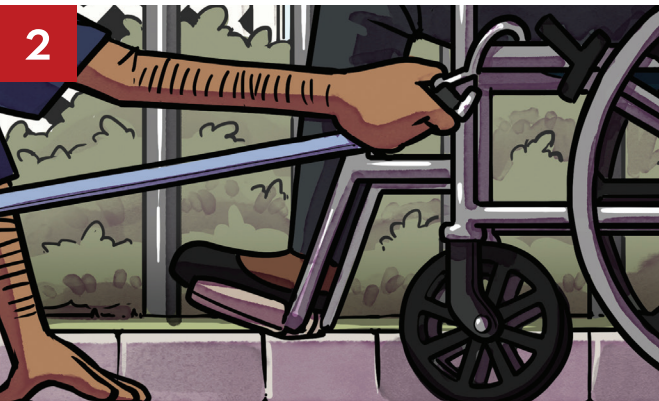


2 Attach hook to securement spot.

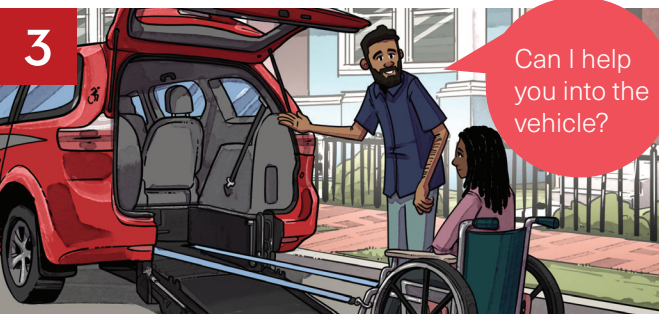


3 Assist passenger into vehicle.

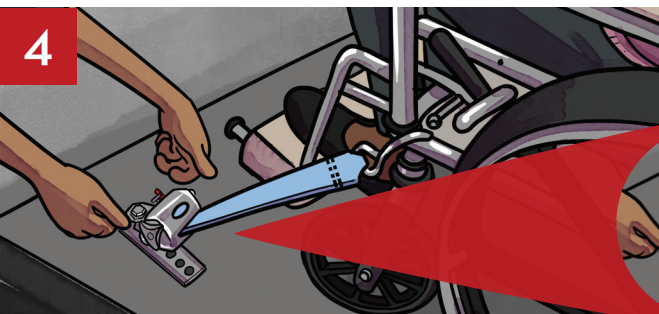


4 Tighten front securements.

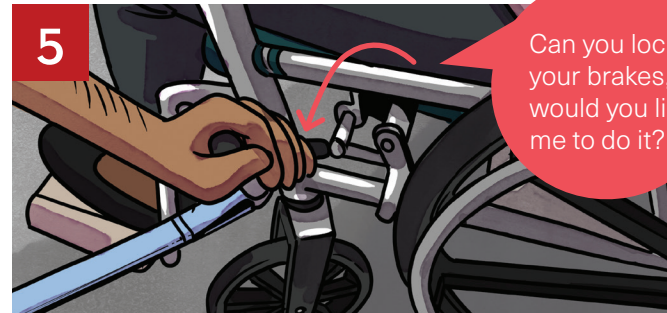

5 Lock the brakes.

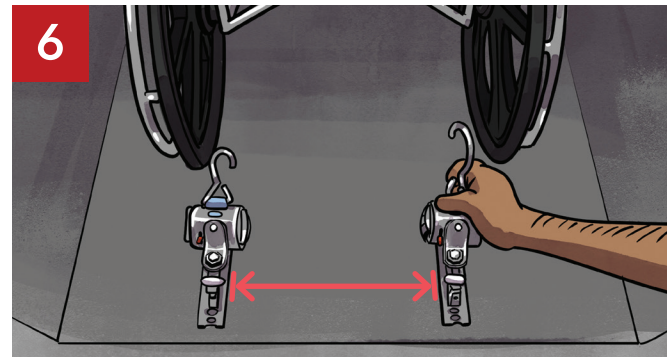


6 Slide back securements into tracks.

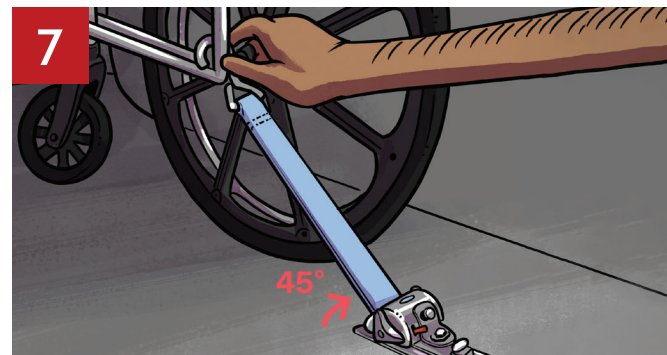

7 Hook straps onto back of wheelchair and tighten.



8 Place lap belt across passenger's lap. The red button should be at right hip.

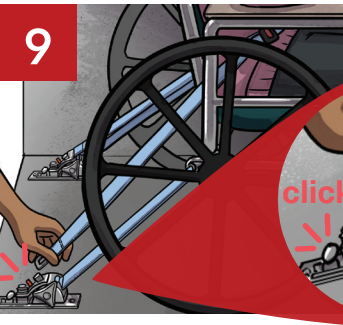


9

Connect lap belt to back securements.



10

Tighten the lap belt.

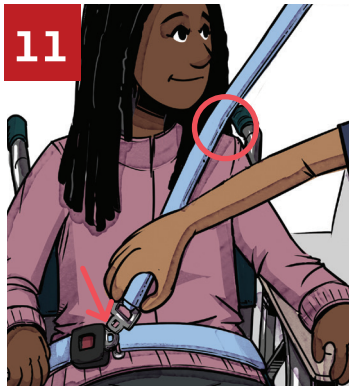

11


Connect the shoulder strap.



12

Tighten the shoulder strap.

Are you comfortable?

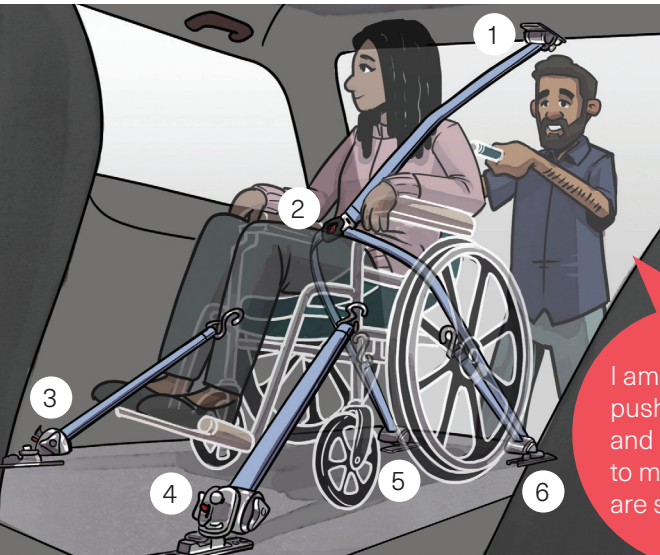

Test the securements by gently pushing and pulling the wheelchair to make sure the securements keep the passenger in place.

I am going to push you back and forth gently to make sure you are secured.

You're Ready to Ride